F

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D

PACKAGE INFORMATION LEAFLET

Aktimed GRID TAPE PLUS is a grid-shaped tape based on polyurethane that has Arnica montana 10⁻⁶ [D6] integrated into the acrylic adhesive. The focus is on the physical and mechanical properties of the adhesive tape.

APPLICATIONS

- Trigger Point and Acupuncture Point Taping *)

INDICATIONS FOR USE

- The product is intended for single use.
- Application is recommended by a trained person such as a doctor, naturopath or physiotherapist.
- Before applying the grid tape degrease the skin.
- The tape should be applied 20-30 minutes before physical activity.
- It may occasionally cause skin irritation. In this case, immediately remove the tape.

SAFETY INFORMATION

- In case of allergy to any of the ingredients do not use the tape or first test on a small area of skin!
- Do not adhere to open wounds or skin injuries!
- In the case of skin reactions such as itching or redness remove the tape immediately!
- Pregnant women may only apply the tape after consultation with the attending physician!
- Keep out of reach of children!

Aktimed GmbH, Theaterstr. 18, 69117 Heidelberg, Germany
www.aktimed-tape.de

E

APPLICATION TIME

The tape lasts up to 6 days on the skin.

INGREDIENTS

80% polyester, 15% polyurethane, 5% acrylate.
Adhesive mass contains Arnica Montana 10⁻⁶ [D6].

TECHNICAL DATA

- **SMALL:** 2.1 cm x 2.7 cm
- **MEDIUM:** 2.8 cm x 3.6 cm
- **LARGE:** 4.4 cm x 5.2 cm

Store in a dry and cool place away from sunlight.
Store between 5° C and 30° C.

CE: conforms to EU Directive 93/42/EEC for medical products.

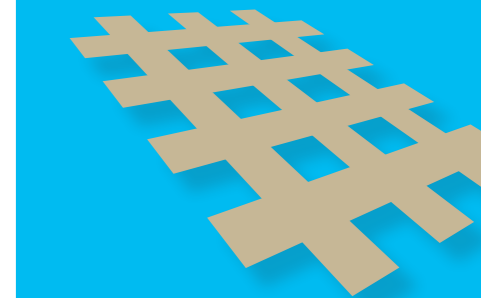
*) Effects have not been scientifically proven

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Aktimed GmbH, Theaterstr. 18, 69117 Heidelberg, Alemania